


FORTUNE ADVANCE ENTERPRISES LIMITED

www.globalsource.com/fadvance.co



Infrared Vibro Mattress

Back lighting LCD command

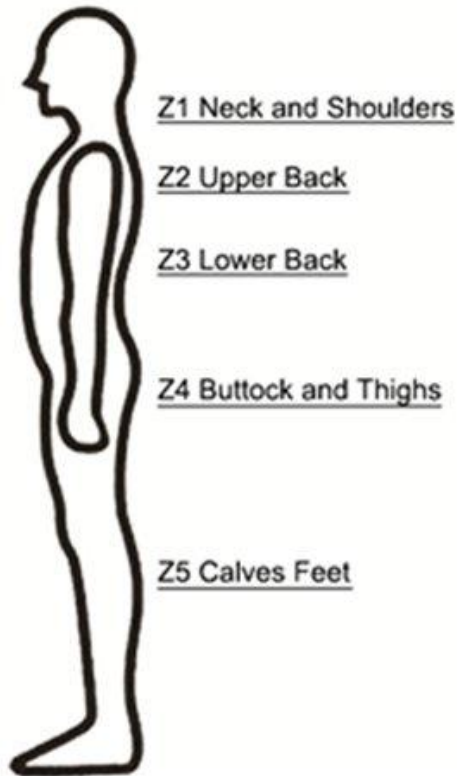

Lavi-660D_

7

- Mattress 195 x 60x 7cm
- LCD command with 24 massager programs
- 7 motors for 5 zones in a mattress
- Each zones works independently
- 7 vibration intensity
- Infrared cushions on pillow, back and foot
- 3 infrared intensity
- Timer
- Massage belt with 2 motors
- Unit price: USD127, 320pcs/ 40



Lavi-660E

9

- Mattress 195 x 60 x 7cm
- LCD command with 24 massager programs
- 9 motors for 5 zones in a mattress
- Each zones works independently
- 7 vibration intensity
- Infrared cushions on pillow, back and foot
- 3 infrared intensity
- Timer
- Massage belt with 2 motors
- Balance cushion
- **UNIT PRICE: USD137, / 259pcs 20'**



Lavi-660F

Mattress 195 x 60 x 7cm

LCD command with 24 massager programs

9 motors for 5 zones in a mattress

Each zones works independently

7 vibration intensity

Infrared cushions on pillow, back and foot

3 infrared intensity

Timer

Massage belt with 2 motors

Balance cushion

UNIT PRICE: USD140, / 259pcs 20'




Balance cushion

Zipped armrest

Lavi-660G_

I3

- Mattress 195 x 73 x 7cm
- LCD command with 24 massager programs
- 13 motors for 5 zones in a mattress
- Each zones works independently
- 7 vibration intensity
- Infrared cushions on pillow, back and foot
- 3 infrared intensity
- Timer
- Massage belt with 2 motors
- Leg balance cushion
- **UNIT PRICE: USD190, /256pcs 40'HQ**

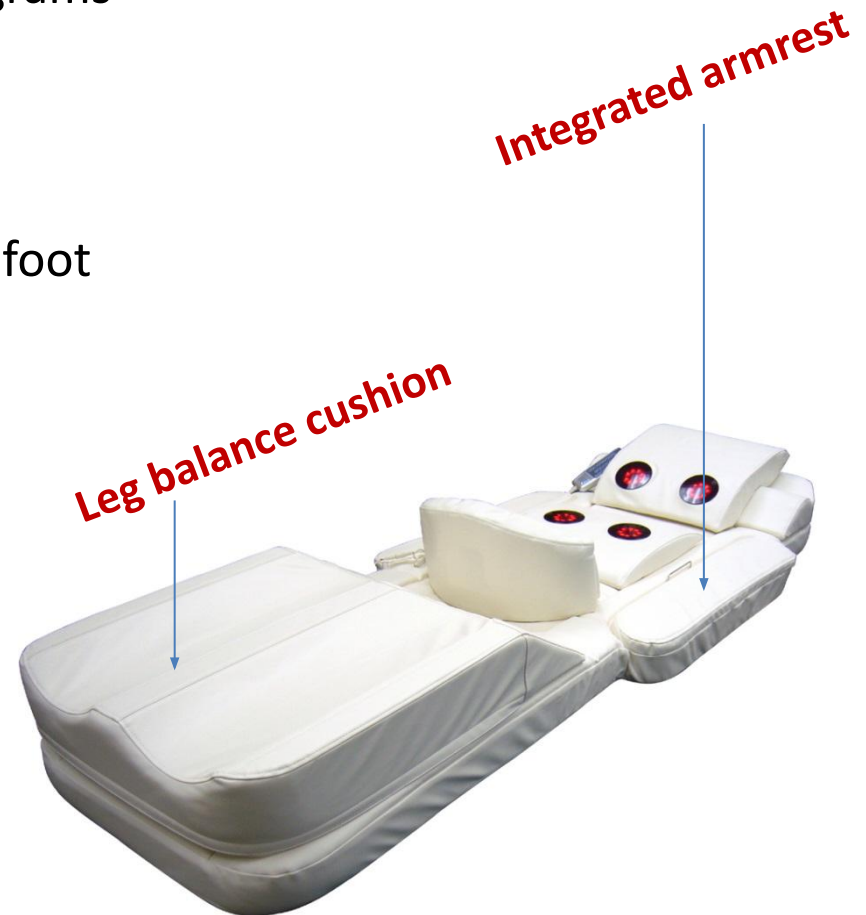

Operation.1

1. First connect all the accessories with the mattress properly.
Second connect the power supply to the mattress then the "ON" LED lights.
The system is ready for work.
2. Select the timer either 5 or 15 or 30 or 60 minutes by function key "T" before
massage operates. If not selected, then 15 minutes will be worked
automatically.
3. Start the massaging by press the function key "P", message program P1
activates. 24 programs are available.
Intensity starts at Z1=1, Z2=3, Z3=3, Z4=3, Z5=3.
-press function key "P" again changes the message program to "P2-P24".
 - 3.1 Adjust intensity for each zone
-press Z1 or Z2 or Z3 or Z4 or Z5 selects the zone first.
-then press "+" "-"adjust the intensity of each message zone individually
-7 levels of intensity are available
-press "-"to 0 turns off each message zone individually
 - 3.2 Adjust Z2-5 zones intensity simultaneously
-press function key "Z2-5" these 4 zones work in same intensity at level 3
simultaneously, it replaces the previous intensity setting.
-then press "+" "-"adjust the intensity of Z2-Z5 zone simultaneously.

Operation.2

4. Infrared heat: (IR1-head cushion, IR2-back cushion, IR3-foot cushion)
 - press function key "IR1", infrared heat 1 start at strong heat intensity.
 - press the key again, infrared intensity changes to medium.
 - press the key again, infrared intensity changes to mild.
 - press the key again, infrared will be turn off.
 - operations are the same for infrared heat "IR2", "IR3".
5. Time
 - press function key "T" to either 5 or 15 or 30 or 60 minutes.
 - while device is in operation, press function key "T" first.
 - then press "+" / "-" function key adjust time in 1 minute.
6. Off
 - press function key "Off" stop the operation at any time.

CERTIFICATE of Conformity



EC Council Directive 2004/108/EC
Electromagnetic Compatibility

Registration No.: AE 50267626 0001

Report No.: 17035737 001

Holder: Fortune Advance Enterprises Limited
Room 2002 Beverly House
93-107 Lockhart Road, Wan Chai
Hong Kong

Product: Massage Equipment
(Massage Mattress)

Identification: Type Designation: Lavi-660
Serial No. : n.a.
Remark: Refer to test report 17035737 001 for details.

Tested acc. to: EN 55014-1:2006+A1+A2
EN 61000-3-2:2006+A1+A2
EN 61000-3-3:2008
EN 55014-2:1997+A1+A2

This certificate of conformity is based on an evaluation of a sample of the above mentioned product. Technical Report and documentation are at the Licence Holder's disposal. This is to certify that the tested sample is in conformity with all provisions of Annex I of Council Directive 2004/108/EC. This certificate does not imply assessment of the production of the product and does not permit the use of a TÜV Rheinland mark of conformity. The holder of the certificate is authorized to use this certificate in connection with the EC declaration of conformity according to the a.m. Directives.

Date 15.10.2013



TÜV Rheinland LGA Products GmbH - Tillystraße 2 - 90431 Nürnberg

CE The CE marking may only be used if all relevant and effective EC Directives are complied with. CE





battery operated
massage cushion



Massage belt
with 2 motors



Massage belt
with 2 motors

Balance cushion